


**Sample Preparation:** 3 pcs. Sorb-Star® in 100 ml water; Restek Nr. 31011,  
24 h / 750 rpm stirred

**Thermal Desorption / VSP4000:** sample: 220 °C; purge time: 20 min;  
purge flow: 20 ml/min;  
trap: -30 °C; trap-desorb: 350 °C


**GC:** HP5, 30m x 0.32 mm x 0,25µm;  
temp. program: 60 °C / 2 min > 20 %/min > 340 °C / 9 min

**Detector:** FID

### Recovery

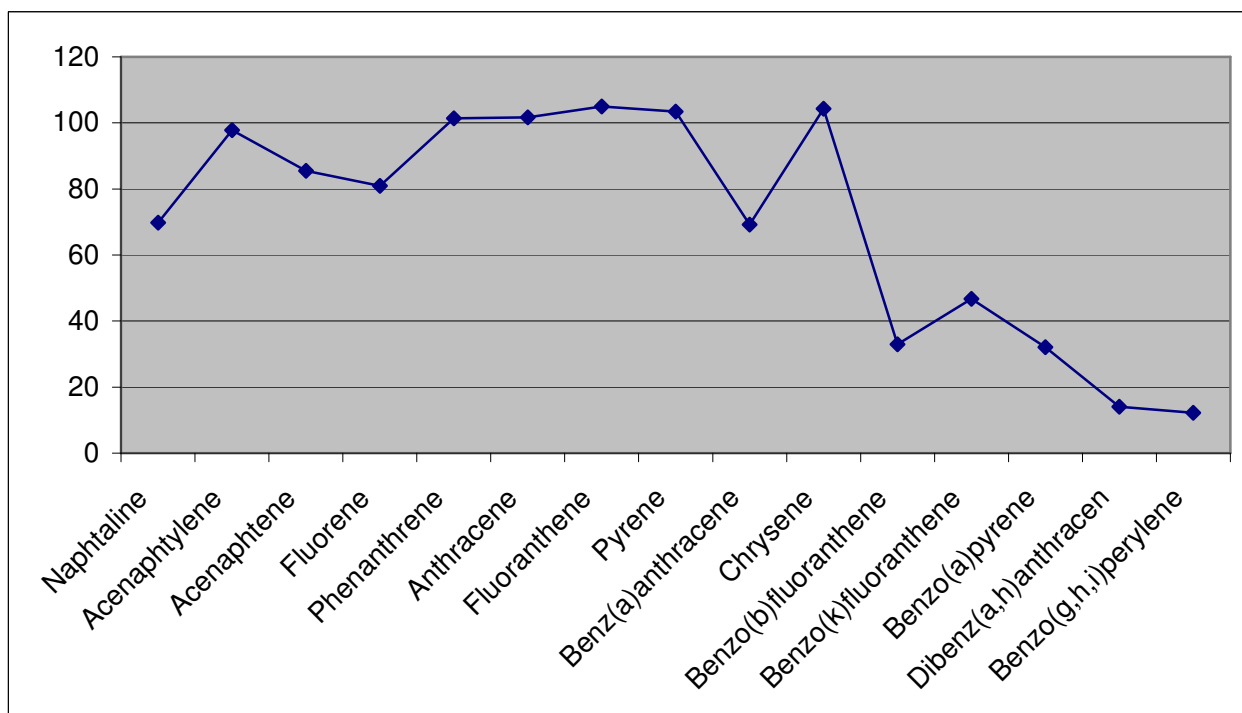

### Chromatogram

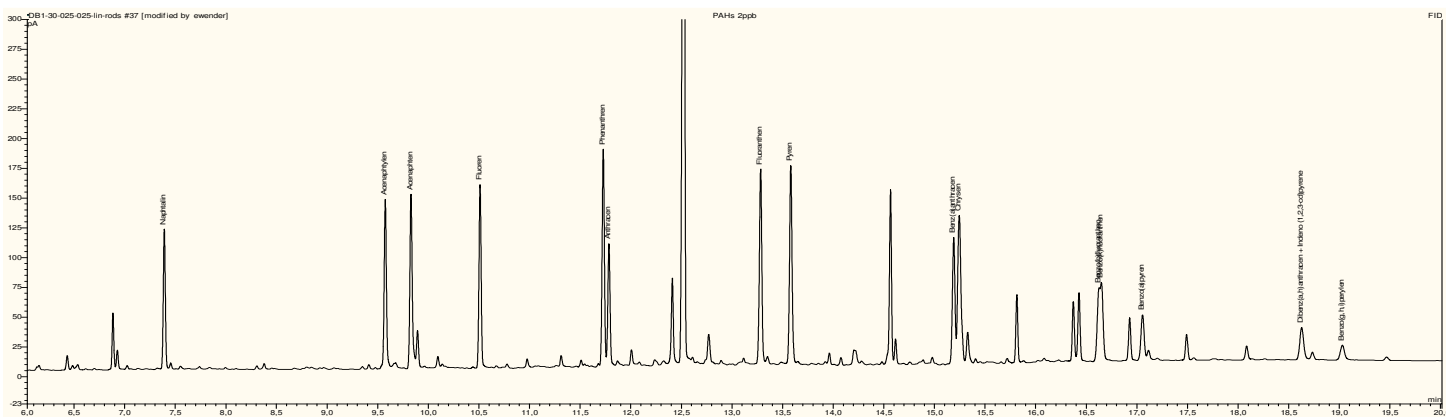

Calibration Curves

